

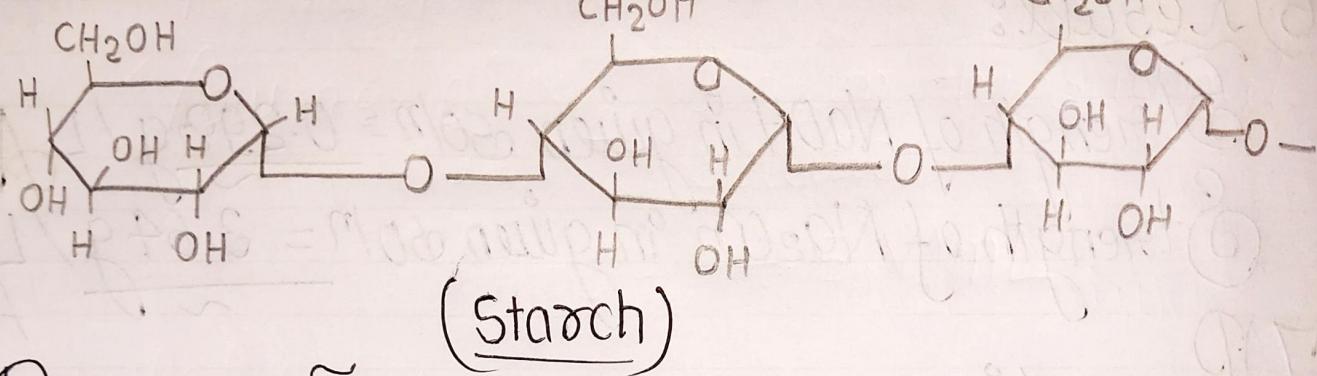
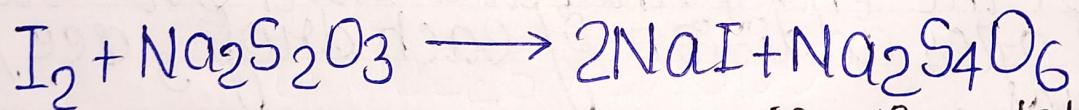
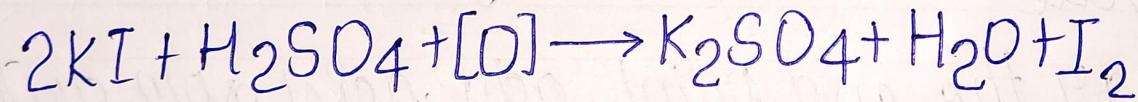
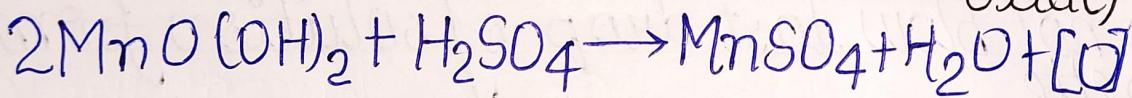


Indicator

Experiment No: A


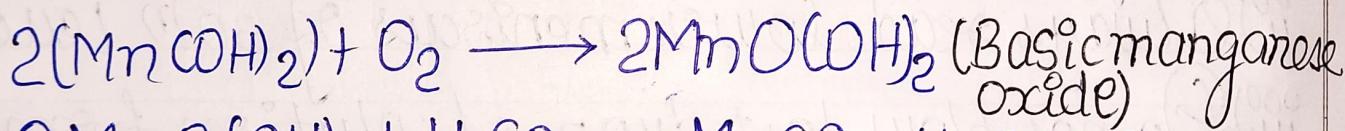


Reactions Involved

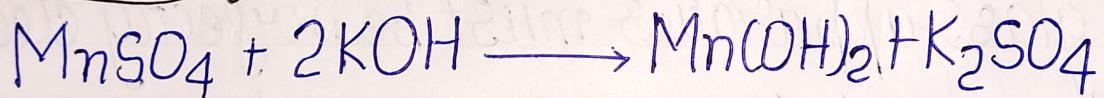

(Sodium tetrathionate)

Tonic Reactions

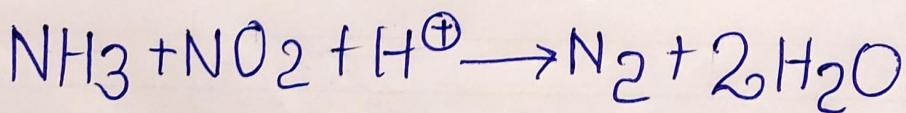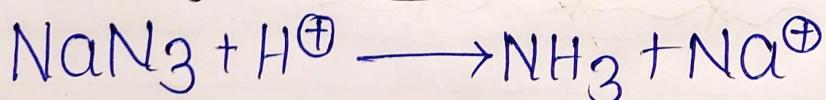

~~bighri~~

* Experiment No: 4 *

- 1) Aim: To determine the concentration of dissolved oxygen from the given water sample using standard N/40 $\text{Na}_2\text{S}_2\text{O}_3$ solution
- 2) Apparatus Required: Burette, Pipette, Measuring flask, Glass Rod.
- 3) Reagents Required: Potassium iodide solution (10%), N/40 Hypo solution ($\text{Na}_2\text{S}_2\text{O}_3$), starch solution (freshly prepared)
- 4) Theory: This experiment is based on the oxidation of potassium iodide by dissolved oxygen. The liberated iodide is against a standard solution. Since, oxygen dissolved in water remains present in molecular state; it is not capable of reacting with potassium Iodide. As oxygen carrier such as manganese hydroxide is produced *in situ* by the actions of potassium hydroxide and manganese sulphate. Starch is used as an Indicator in Iodometric titration. It is a visible indicator in titration process because it turns deep dark blue when iodide is present in a solution when starch is added.

heated in water, decomposition occurs and beta-
amylase is produced Beta-amylase with
iodine, resulting in a dark blue coloured
change.

Iodide modification is used to remove
substances like nitrates, sulphide, etc.
which liberate iodine from potassium
iodide to dissolve oxygen.

For eg: the interference of nitrite is overcome
by adding sodium iodide (NaI) and
 H_2SO_4 .

5.) Procedure:

- 50mL of water sample was taken in Iodometric flask, avoiding as far as possible contact with air.
- Immediately 1ml of manganese sulphate solution was added from the burette and add 1mL of alkaline iodide was added from other burette
- The stopper was inserted and solution was shaken several times

The precipitate was settled half way & was mixed again. This process was repeated at least

Observations :-

Volume of water used = 50 ml

Volume of (N/40) sodium thiosulphate solution

$$V_{\text{Hypo.}} = 44.5 (\text{V}) \text{ mL}$$

Calculations :-

$$N_{\text{O}_2} \times N_{\text{water}} = N_{\text{Hypo}} \times N_{\text{Hypo}}$$

$$N_{\text{O}_2} = \left(\frac{1}{40} \right) \times \left(\frac{44.5}{50} \right) = \underline{\underline{0.02225}}$$

$$\begin{aligned}\text{Strength} &= N_{\text{O}_2} \times 8 \text{ gm of O}_2 / \text{litre} \\ &= N_{\text{O}_2} \times 8 \times 1000 \text{ ppm} \\ &= 0.02225 \times 8 \times 1000 \\ &= \underline{\underline{178 \text{ ppm}}}\end{aligned}$$

~~178 ppm~~

for 3 times.

- 6) The yellow solution was allowed to stand for 5 min.
- 7) 50 mL of solution was withdrawn. It was then titrated against the standard N/40 sodium thiosulphate solution till the solution became faint yellow.
- 8) At this stage few drops of starch were added. The colour turned to deep blue information of starch iodide.
- 9) Titration was continued till the solution becomes colourless.

6) Result: Amount of Dissolved Oxygen from the water using sample N/40 $\text{Na}_2\text{S}_2\text{O}_3$ solution standard

$$= 178 \text{ ppm.}$$

7) Precautions:

- (i) The solution in the iodide flask should be shaken carefully while keeping the thumb on the lid.

- (iii) The solution be well shaken before each aliquot is withdrawn for titration
- (iv) The solution in the iodine flask should not be sucked into the pipette with the mouth

~~bjohari~~