

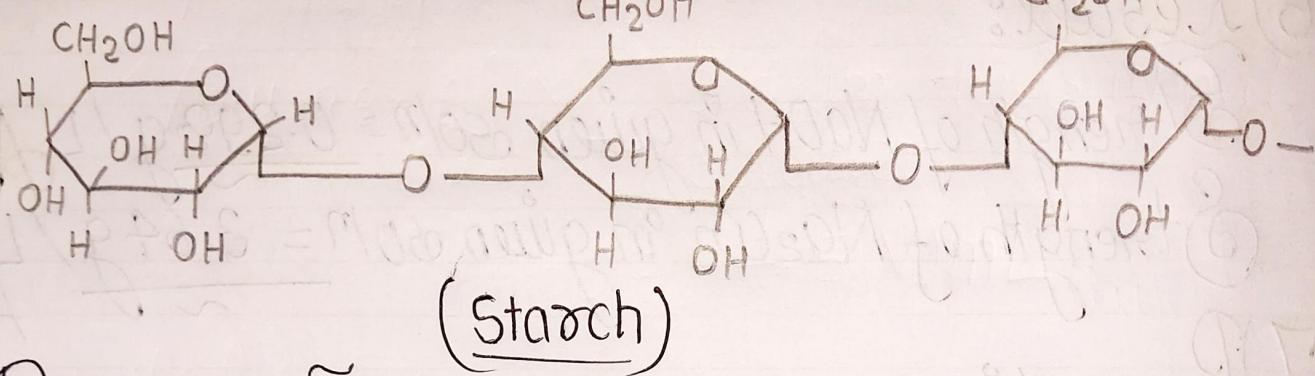
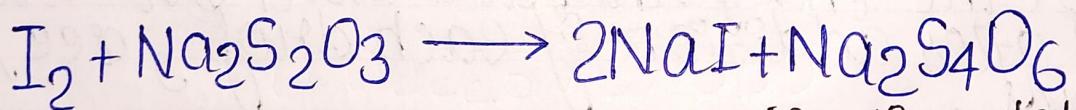
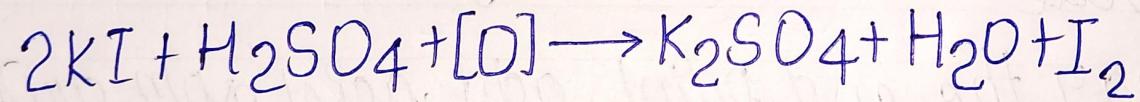
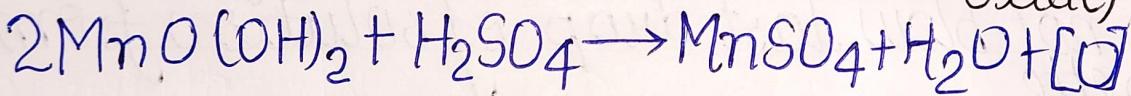


## Indicator

## Experiment No: A

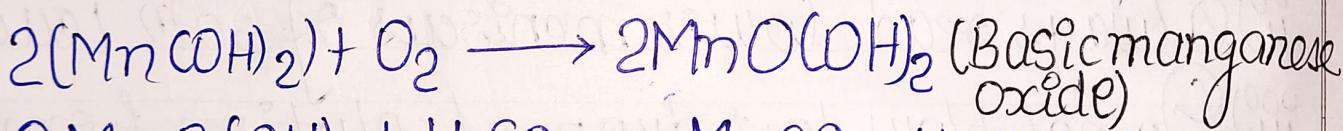


## Reactions Involved

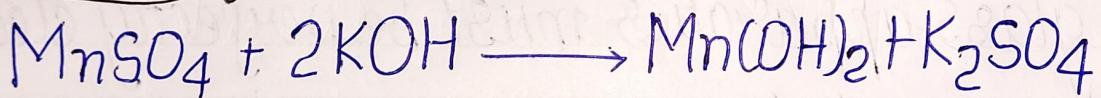

(Sodium tetrathionate)

## Tonic Reactions

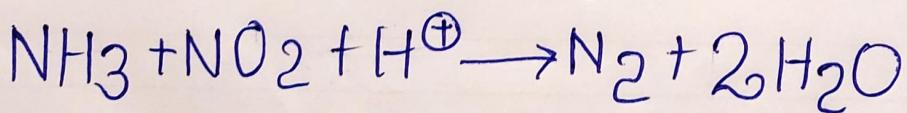

~~bighri~~

# \* Experiment No: 4 \*

- 1) Aim: To determine the concentration of dissolved oxygen from the given water sample using standard N/40  $\text{Na}_2\text{S}_2\text{O}_3$  solution
- 2) Apparatus Required: Burette, Pipette, Measuring flask, Glass Rod.
- 3) Reagents Required: Potassium iodide solution (10%), N/40 Hypo solution ( $\text{Na}_2\text{S}_2\text{O}_3$ ), starch solution (freshly prepared)

4) Theory: This experiment is based on the oxidation of potassium iodide by dissolved oxygen. The liberated iodide is against a standard solution. Since, oxygen dissolved in water remains present in molecular state; it is not capable of reacting with potassium Iodide. As oxygen carrier such as manganese hydroxide is produced *in situ* by the actions of potassium hydroxide and manganese sulphate. Starch is used as an Indicator in Iodometric titration. It is a visible indicator in titration process because it turns deep dark blue when iodide is present in a solution when starch is added.

heated in water, decomposition occurs and beta-  
amylase is produced Beta-amylase with  
iodine, resulting in a dark blue coloured  
change.

Iodide modification is used to remove  
substances like nitrates, sulphide, etc.  
which liberate iodine from potassium  
iodide to dissolve oxygen.

For eg: the interference of nitrite is overcome  
by adding sodium iodide ( $\text{NaI}$ ) and  
 $\text{H}_2\text{SO}_4$ .

### 5.) Procedure:

- 50mL of water sample was taken in Iodometric flask, avoiding as far as possible contact with air.
- Immediately 1ml of manganese sulphate solution was added from the burette and add 1mL of alkaline iodide was added from other burette
- The stopper was inserted and solution was shaken several times

The precipitate was settled half way & was mixed again. This process was repeated at least

## Observations

Volume of water used = 25 mL

Volume of N/40 thiosulphate solution  
= 33.4 mL

## Calculations

25 mL of water sample  $\approx$  V mL of N/40 Na<sub>2</sub>S<sub>2</sub>O<sub>3</sub> solution

$$N_{H_2O} \times V_{H_2O} = N_{Hypo} \times V_{Hypo}$$

$$N_{H_2O} = \frac{N_{Hypo} \times V_{Hypo}}{V_{H_2O}}$$

$$= \left( \frac{1}{40} \right) \times \frac{33.4}{25}$$
$$= 0.0334 N$$

$$\text{Strength (in ppm)} = N \times E_{eq} \text{ l/l} \times 1000$$
$$= 0.0334 \times 9 \times 1000$$

$$= 300.6 \text{ ppm}$$

$\therefore$  The amount of Dissolved

Oxygen = 300.6 ppm ~~higher~~

for 3 times.

- ⇒ The yellow solution was allowed to stand for 5 min.
- ⇒ 50 mL of solution was withdrawn. It was then titrated against the standard N/10 sodium thiosulphate solution till the solution became faint yellow.
- ⇒ At this stage few drops of starch were added. The colour turned to deep blue information of starch iodide.
- ⇒ Titration was continued till the solution becomes colourless.

6) Result: Amount of Dissolved Oxygen from the water using sample N/10  $\text{Na}_2\text{S}_2\text{O}_3$  solution standard  
= ~~178~~ ~~300.6 ppm~~ Ans

7) Precautions:

- (i) The solution in the iodide flask should be shaken carefully while keeping the thumb on the lid.

- (iii) The solution be well shaken before each aliquot is withdrawn for titration
- (iv) The solution in the iodine flask should not be sucked into the pipette with the mouth

~~bjohari~~